

Action 5

THE PERFECT FIT FOR YOUR LIFESTYLE



The Invacare Action⁵ has been designed for the more active user looking for an affordable chair with good drive performance and the ability to fold easily for handling, storage and transportation purposes. With its two patents pending and lightweight construction, the new Action⁵ gives the user a chair they can rely on.



Yes, you can.®

Making Life's
Experiences Possible

Action 5



Designed specifically for users with an active lifestyle, the Invacare Action 5 Folding is a highly adjustable, stable chair with a unique and innovative folding mechanism that allows the backrest to fit neatly on top of the seat.


The Action 5 features a light weight frame with step-less adjustability, meaning very fine adjustments can be made. In addition many of the chair's features can be adjusted while the user remains seated. Swing away and removable armrests allow the Action 5 to be reconfigured for a variety of different situations.

The Action 5 is available in a wide range of attractive colours as well as three accent colours for smaller parts like the castors and rear wheel attachments.

Features



Tension and depth adjustments can be made via easy to use straps.

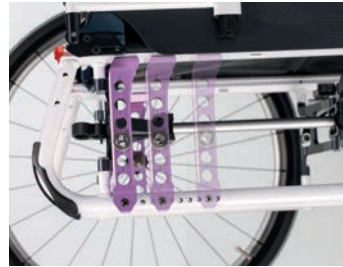

100mm centre of gravity adjustments can be made in 25mm increments which can even be fine-tuned to 12mm.



Standard adjustments of 0-2° available with the option of 3 & 4° adjustments if required.



A wide range of options and accessories are available to give your Action 5 the personal touch.

Tech Specs



Seat width

30.5 / 33 / 35.5 /
38 / 40.5 / 43 /
45.5 / 48 / 50.5 /
53 cm



Seat depth

35 / 37.5 / 40 /
42.5 / 45 / 47.5 /
500 cm



Seat height

front: 36 - 53.5cm
rear: 36 - 50cm
(increments of
2.5cm)



Total width

SW +16 cm To
SW+ 25 cm



Total width folded

32 -38 cm



Total length

84.5 -120.5 cm



Max. user weight
100 & 130 kg



Transport weight
w/out rear wheels
8.3 kg



Total weight
unloaded
from 12.3kg



Turning radius
105 -138 cm



Max. safe slope
7°



Crash Test
Approved
Fulfills the
requirements
of ISO 7176-19

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 05/2024



Yes, you can.®