

# Küschall®



**| UNLIMIT YOUR WORLD**

KSL | K-SERIES | COMPACT | CHAMPION



Yes, you can.®

**Making Life's  
Experiences Possible**





# ***UNLIMIT YOUR WORLD***

We endeavor to offer superior solutions for you. We believe anything is possible when you have the power to achieve it. We want you to live an active life and have the confidence and drive to not give up in the face of any adversity.

We are determined to offer the best solutions. We will Unlimit your world.



# **HYDROFORMING**

Hydroforming is a process which allows metal to be formed into complex shapes. This enables the creation of more ergonomically designed pieces, which give our chairs a unique look and provide improvements in user comfort through better ergonomics.

Hydroforming also allows variances in wall thickness to optimize the ratio of weight versus strength and rigidity. By using hydroforming technology, the Kuschall frames are not only stronger but also allow more effortless propulsion. Tests show that for the same effort generated on the frame, deformation on a hydroformed frames is 60% lower than the rounded tube design and equally comparable to the performance of carbon fibre. Put simply; Hydroforming helps create a wheelchair that is stronger than before and effortless to propel.



# THE KSL®

| RIGID



Elegance is at the core of the KSL®, get your full Kuschall experience with unrivalled quality and minimalistic design. Designed for the more experienced user, the made-to-measure KSL® is fully customisable giving you the freedom to express your own personality. The Super Light KSL® is the essence of pure beauty. With the fully welded hydroformed frame, the rigidity provides a graceful effortless glide anywhere you go.

- Elegant
- Made-to-measure
- Super Light



Yes, you can.®



SHIMANO  
KINOSHITA





1



2



3



4



1. Simply the lightest  
Made-to-measure concept for the lightest construction. Kuschall Super Light brings all of that plus a powerful geometry that flows seamlessly into the chair's silhouette.
2. Integrated design  
From the front frame, to the back brace and the fork support, the new KSL® speaks the same design language. Your new KSL® has been powerfully and clearly designed, uniting traditional and modern.
3. Fix or foldable backrest  
While being designed to enhance stiffness and strength, KSL® can also be ordered with the foldable backrest. Never settle for compromise and take it anywhere you go. The KSL® is perfect for transport and storage.
4. Re-designed footrest  
Details always make the difference. We have re-designed the footrest, a fundamental element for your seating posture, following the powerful geometry and dynamic appeal of the overall look.

# K-SERIES®

| RIGID



The iconic modular K-Series is the perfect product for you. The unique versatile seat module can be effortlessly modified independently from the frame and is fully adjustable to be able to fit you perfectly. Whether you require a more stable position or are more advanced, the K-Series can be adapted to suit your lifestyle. The dynamic K-Series offers you a superb driving performance that will keep up with you on the move.

- Iconic
- Modular
- Dynamic



Yes, you can.®



1



2



3



4



# K

1. Timeless innovation  
This chair is the essence of Küschall and what the brand stands for, including all the K-series know how that users want. A wheelchair that combines versatility, performance and design. If there's a product that represents everything we believe in, K-Series is it.
2. Seat module  
If you are looking for another chair like this, you won't find it. The seat module remains independent from the frame and allows you to change your seating set-up anytime, anywhere.
3. Ergonomic backrest brace  
More grip and less lifting effort thanks to the new sleek hydroformed backrest brace. A cleverly designed wheelchair, engineered to adapt to the changing elements in your daily life.
4. Quick and easy folding  
The lockdown folding backrest is height and angle adjustable. Easy folding with the perfect open frame design ensures effortless car transfers. Keep independent and on the go!

# CHAMPION®

| FOLDABLE




This hybrid product is foldable with the rigidity of a fixed frame giving you an incredible driving performance. The Champion is the perfect companion for your adventures. The chair becomes compact for travel using the unique folding mechanism and with the SK (Swiss knife-folding) option the Champion can be folded even smaller. With the innovative lightweight hydroformed frame the Champion is ideal for your exciting life on the move.

- Unique
- Hybrid
- Super Compact



Yes, you can.®

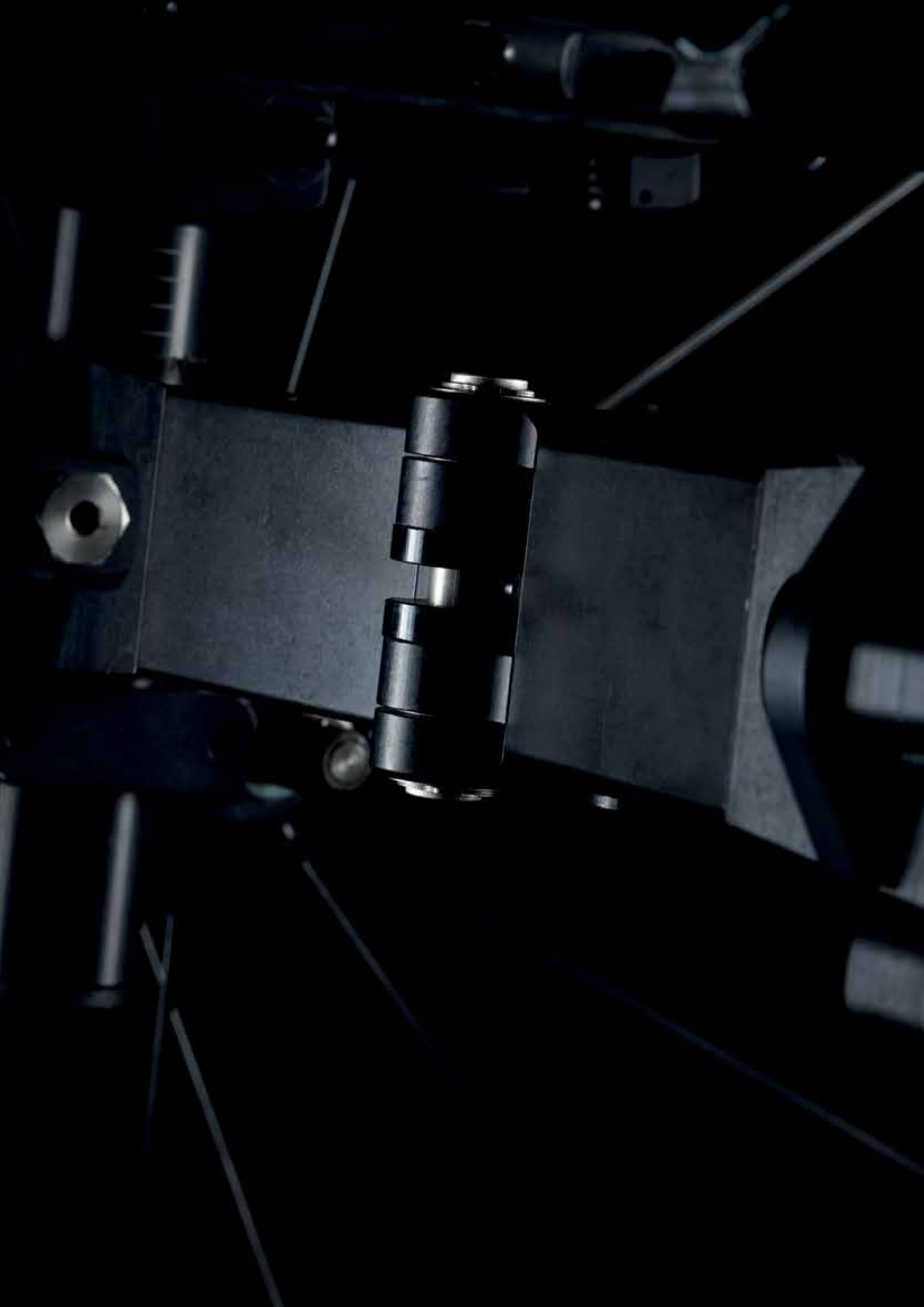





1

2



3

4

**K**

1. Unique folding mechanism  
Unique, because nothing as perfect can be found on the market. The locking mechanism is at the heart of the Champion. All the qualities of a folding wheelchair and the necessary stiffness to achieve the driving performance of a rigid frame. Hydroforming gives the final touch, because a Champion always strives for greatness.
2. Sleek redesign of the one-piece flip up footplate  
The new footrest flawlessly adapts to the new design language of the portfolio. Clear-cut shapes and contours to give the best geometry to the whole wheelchair, keeping its technical benefits unaltered.
3. Small to the size of a backpack  
If traveling is your dream, Champion Swiss Knife option will make it come through. The SK option features a folding front frame, easy to close, that integrates perfectly with the lines obtained with hydroforming. A new generation, a new beginning, for an even better Champion.
4. New light active brake  
Designed with modern breakthrough lines and lighter by almost one third compared to other models, the new brake saves weight and gives you perfect control during activation. Precision even in the smallest details.

# COMPACT®

# | FOLDABLE




The Compact is the absolute standard for reliability, high-class quality and outstanding functionality. Küschall Compact can easily be adapted to suit your needs and lifestyle, whether you need a higher level of configurability or an active product, the Compact is right for you! This traditional style, cross folding wheelchair can be effortlessly folded. The lightweight frame is robust giving you a smooth driving performance.

- Reliable
- Outstanding Functionality
- Enhanced Comfort



Yes, you can.®




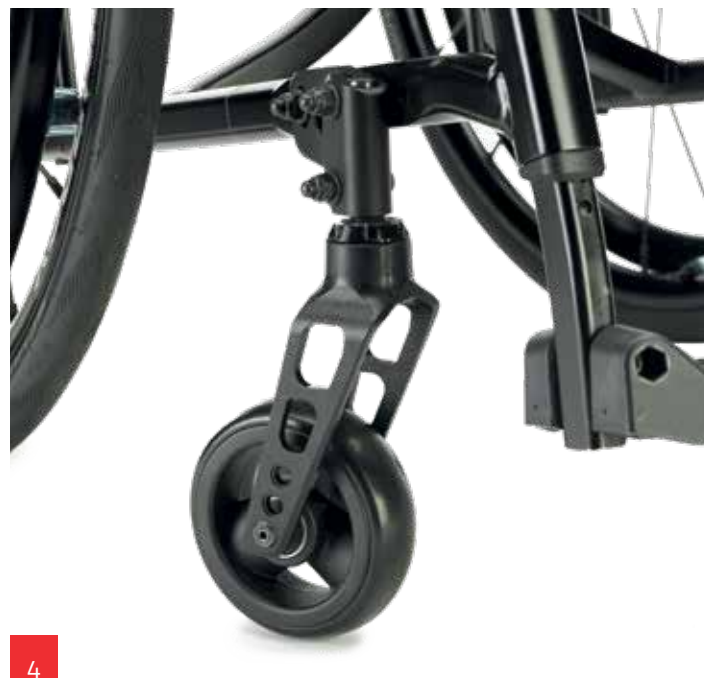

1. Available with fixed or swing-away leg rests.  
We have wisely merged versatility and performance. With the new Küschall Compact you can choose the fixed or swing-away front frame, the result is an all-around product with the stance of a super active wheelchair.
2. More stiffness with smaller tubes  
Hydroforming allows a significant reduction of the weight of the frame, guaranteeing less effort in propelling, easy front transfers and a compact storage. The combination of design and built-in precision results in an ultra-light wheelchair that delivers on performance.
3. Ergonomic design  
Thanks to the new shapes in the knee area you will benefit from maximum grip for easy transfer, repositioning or loading in the car. Ergonomics meets today's hydroforming possibilities, giving life to a design that unites tradition and modernity.
4. Keep on track  
Your needs drive our innovations. We have designed this intuitive castor housing adjustment, which is excellent in keeping the direction steady'.

# CUSTOMISATION

## 1 CHOOSE THE FRAME COLOR



Pure White



Anthracite Grey



Space Black



Azure Blue



Carmine Red

## 2 CHOOSE THE DECAL STYLE



## 3 CHOOSE THE DECAL COLOR



Black



Red



Light Grey



Gold



Red



Blue



Dark Grey



Black

## KÜSCHALL VISUALISER




Connect to [www.kuschall.com/international-en/visualizer/](http://www.kuschall.com/international-en/visualizer/) and create something unique with the Kuschall Visualiser. Design your wheelchair around your style and break the boundaries!

# ACCESSORIES

Do not set limits to your desires, we have what you need and much more.

Go to [www.kuschall.com/accessories](http://www.kuschall.com/accessories) to explore the entire range of accessories for your chair.

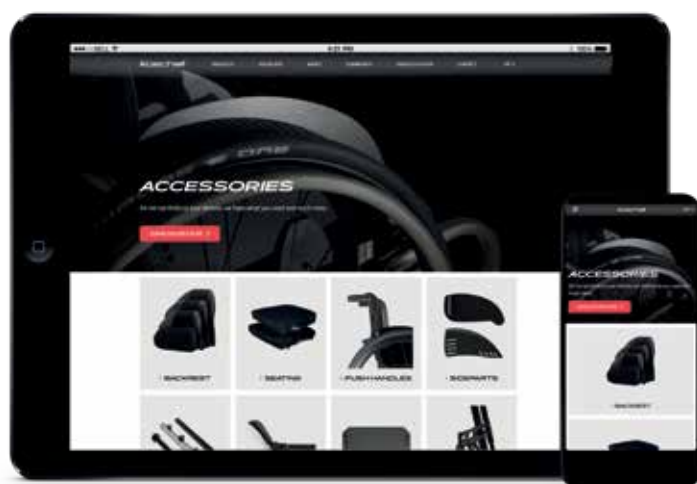

# TECH SPECS

## RIGID

	<b>THE KSL®</b>	<b>K-SERIES</b>
▪ Frame angle	75° / 90°	75° / 90°
▪ Seat width	320 - 500 mm	320 - 500 mm
▪ Seat depth	350 - 525 mm	350 - 525 mm
▪ Total width	Seat width + 170 mm	Seat width + 170 mm
▪ Total length	75°: approx. 800 mm 90°: approx. 730 mm	75°: approx. 830 - 1055 mm 90°: approx. 760 - 985 mm
▪ Total weight unloaded	approx. 6.8 kg	approx. 7.8 kg
▪ Transport weight (w/out rear wheels)	approx. 4.3 kg	approx. 5.1 kg
▪ Max. user weight	130 kg	130 kg

## FOLDABLE

	<b>CHAMPION</b>	<b>COMPACT</b>
▪ Frame angle	75° / 90°	70° / 80° (swing away, fixed and tapered)
▪ Seat width	360 - 480 mm	280 - 500 mm
▪ Seat depth	340 - 480 mm	320 - 500 mm
▪ Total width	Seat width + 160 -240 mm	Seat width + 160 mm
▪ Total length	75°: approx. 930 mm 90°: approx. 830 mm	70°: approx. 870 – 1220 mm 80°: approx. 850 – 1200 mm
▪ Total weight unloaded	approx. 9.2 kg	approx. 10.5 kg (swing away) approx. 9.6 kg (fixed frame)
▪ Transport weight (w/out rear wheels)	approx. 6.7 kg	approx. 7.8 kg (swing away) approx. 6.7 kg (fixed frame)
▪ Max. user weight	120 kg	130 kg

All measurements and weights indicated are based on seat width 400 mm and seat depth 400 mm in the lightest and most compact configuration.

Weight and measurements can change with the different configurations. For further information please check the user manual on our website.

**K** Kuschall is a registered trademark.  
Every effort has been made to ensure that the contents of this publication are fully up-to-date at the time of going to print. As part of its ongoing improvement of products, Kuschall AG reserves the right to modify existing models at any time. Any use of this publication or of parts thereof as well as any reproduction of images must have the written consent of Kuschall AG.  
2019 Copyright © 11/2019, Kuschall AG, Switzerland. All rights reserved.



**Yes, you can.®**



FOR FURTHER INFORMATION:

**New Zealand:** Freephone 0800 468 222, email [sales@invacare.co.nz](mailto:sales@invacare.co.nz) or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:** Freephone 1800 460 460, email [orders@invacare.com.au](mailto:orders@invacare.com.au) or visit [www.invacare.com.au](http://www.invacare.com.au)

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 07/2024



**Yes, you can.®**