

Softform Premier

DURABILITY AND PROTECTION



Equipped with a suite of features that ensures durability, sustained high degree of protection against tissue damage for the patient and an ingeniously designed insert pad, the Softform Premier mattress really does provide exceptional pressure reduction and comfort.

High-Very High Risk



Yes, you can.®

Making Life's
Experiences Possible



Softform Premier


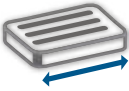
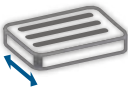
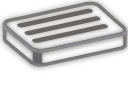
Unique Features

Benefits

Deep contoured insert pad for enhanced pressure care area	Maximises body contact area and optimises pressure reduction
Cover fabric conforms to deep contoured foam castellations on the upper surface of the insert pad when mattress is profiled	Minimizes rucking of cover fabric on top surface, reducing potential for tissue damage
Suitable for High Risk patients / Very High Risk when used in conjunction with a profiling bed	Can contribute significantly to pressure ulcer prevention
No need to turn over	Reduces manual handling issues
Single piece foam 'U' core	Facilitates patient transfer by providing effective side-wall support
Double-coated poly-urethane sides	Reduces the risk of strike through, hydrolysis or delamination
Multi-way stretch polyurethane waterproof cover	Provides excellent recovery and durability and helps reduce shear and friction forces
High frequency welded seams and zip	Reduces risk of fluid ingress
Full length waterfall flap conceals the zip	Further risk reduction of fluid ingress (through zip)
Toughened PU coated base	Prolongs longevity of the mattress
Water penetration resistance	35 kPa minimum (BS 3424: part 26). 35 kPa minimum after 20,000 cycles Martindale Abrasion Test EN ISO 5470-2

Tech Specs

Codes

					WARRANTY
	User weight	Length	Width	Depth	Warranty
Single	247.6kg	1970mm	880mm	152mm	4yr cover 8yr foam
King Single	247.6kg	1970mm	1070mm	152mm	

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 07/2024



Yes, you can.®