

AQUATEC®

# Ocean Ergo

adjustments and options



Simple adaptations that make a difference

Adjustable, ergonomic and versatile



Yes, you can.®

# Adjustments and options

Our Ocean Ergo range offers an advanced solution in shower chairs. Each chair is designed to offer high levels of adjustments and configurability as standard. For further customisation, our latest range of options provides solutions for extra comfort, support and safety.



### Footplates

- > Angle adjustable, flip-up and textured footplates with heel loops
- > No tools required

### Footrests

- > Height adjustable, swing away and detachable footrests
- > No tools required

### Armrests

- > Ergonomic, swing away armrests can be adjusted to different heights

### Seat height

- > Seat to ground height can be adjusted between 440-600 mm\* or 475-600 mm\*\*
- \*Ocean Ergo, Ergo XL and 24" Ergo models
- \*\*Ocean 24" Ergo XL and VIP models
- > No tools required

# Options

## Seating and backrests



**Soft seat**  
Redesigned to mirror the contours of the new ergonomic seat plate, the softer feel gives users with sensory or skin protection issues added relief. Smaller aperture also available for shorter or thinner users. Easy to fit and clean.



**Soft seat insert**  
Can be used with the soft seat when the commode function is not required. Available in two versions for standard soft seat and small soft seat. Standard version fits also directly on the seat plate.



**Ergonomic soft seat**  
Manufactured from viscoelastic foam, it offers a much softer feel for those with skin integrity or sensory issues. The foam creates greater depth of immersion promoting increased pelvic stability and positioning.



**Variable soft seat**  
Hygiene cut out can be positioned in all four directions, improving access for carers assisting with personal hygiene.



**Universal soft seat**  
Comfortable, closed seat surface for when commode function is not in use.



**Solid backrest**  
The solid back support offers an alternative to the standard tension adjustable backrest, when firmer support is required to achieve a more upright posture. Easy to clean.



**Armrest pads**  
New design offers improved support and comfort for the user (not for XL armrests).



**XL back**  
Increases the width between the armrests by 80 mm. Can be retrofitted to standard Ocean Ergo models.



**Wider XL armrests**  
Increases the distance between the armrests by 40 mm on each side.



**Soft backrest cushion**  
For increased comfort. Can also be used to reduce the overall seat width and depth, helping to position smaller users.



**Comfort backrest**  
An alternative to the standard backrest with a softer feel material for increased comfort.



**Ocean VIP Ergo backrest frame**  
For retrofitting Ocean Ergo or Ocean 24" Ergo in order to use a headrest or neckrest (headrest, armrests, backrest cover not included).

## Safety



**Safety bar**  
For additional safety and reassurance.



**Safety lock for armrest**  
Locks armrest into the vertical or horizontal position.



**Front and rear anti-tippers**  
For additional safety plus front and rear stability.



**Tilt-stop**  
An aid for those who are hoisted in/out of the chair, the tilt-stop screws ensure the frame remains horizontal.

## Positioning



**Safety belts**  
Broad pad manufactured from soft material with single pull buckle. Can be used across the pelvis or in addition across the chest for positioning and safety. Easy to attach. Machine washable.



**Lateral supports**  
Adjustable in height via three possible positions on the backrest cane. Broad lateral pads can be repositioned for depth adjustment, comfort and effective pressure redistribution.



**Ergonomic headrest holder**  
Offers multidirectional adjustment to correct or accommodate for more complex head/neck positioning needs. Pad not included.




**Ergonomic headrest/neckrest**  
Offers exceptional management of the head position. Neckrest shown.

## Toileting



**Calf and amputee support with foot support**  
Offers both height and angle adjustment with a contoured pad for added positioning and comfort. Including a foot support for those who may need to elevate their lower limb(s) for medical reasons or comfort.

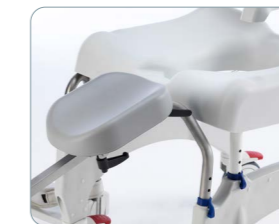

**Calf and amputee support**  
Excluding footrest.



**Splash guard**  
Clips on easily for use with the standard plastic seat, standard soft seat and small soft seat.



**Toilet pan**  
Ergonomic handle allows the pan to be easily removed and replaced. Features an integrated lid holder, so the lid is not misplaced during use.

# Customise each of our Ocean Ergo chairs

Our Ocean Ergo range of shower chairs can all be customised to suit a wide number of needs. By making a simple adjustment or by adding one of our easy-to-fit options, each of our modular shower chairs are more versatile than ever before.



### Ocean Ergo

An entry level attendant propelled model



### Ocean Ergo XL

An entry level attendant propelled model with a heavier maximum user weight



### Ocean 24" Ergo

An entry level self-propelling model



### Ocean 24" Ergo XL

An entry level self-propelling model with a heavier maximum user weight



### Ocean VIP Ergo

A tilt-in-space model for additional seating and positioning support



### Ocean Dual VIP Ergo

A tilt-in-space model with recline for additional seating and positioning support

## Compatibility

✓ Compatible  
 x Not compatible  
 + Included

Most options are cross compatible with all Ocean Ergo models.

	Product code	Ocean Ergo	Ocean 24" Ergo	Ocean Ergo XL	Ocean 24" Ergo XL	Ocean VIP Ergo	Ocean Dual VIP Ergo
Standard soft seat	AP1603142	✓	✓	✓	✓	✓	✓
Soft seat insert	AP1603143	✓	✓	✓	✓	✓	✓
Small soft seat	AP1603145	✓	✓	✓	✓	✓	✓
Small soft seat insert	AP1603152	✓	✓	✓	✓	✓	✓
Splash guard	AP1603155	✓	✓	✓	✓	✓	✓
Universal soft seat	AP1603154	✓	✓	✓	✓	✓	✓
Ergonomic special soft seat	AP60121311	✓	✓	✓	✓	✓	✓
Variable soft seat	1544469	✓	✓	✓	✓	✓ <sup>1</sup>	✓ <sup>1</sup>
Safety bar 2.0 (only with armrests 2.0)	AP60123732	✓	✓	x	x	✓	✓
Safety bar 2.0 XL (only with armrests 2.0)	AP60123771	x	x	✓	✓	x	x
Safety bar XL (only with wider XL armrests)	1470078	✓	✓	x	x	✓	✓
Chest belt	1470081	✓	✓	✓	✓	✓	✓
Pelvic belt	1470082	✓	✓	✓	✓	✓	✓
XL back	AP1604015	✓	✓	+	+	x	x
Solid backrest	1575299	✓	✓	x	x	✓	✓
Soft backrest cushion	1558769	✓	✓	✓	✓	✓	✓
Armrest pad (1 pad; only for new armrests 2.0)	AP1643748	✓	✓	✓	✓	✓	✓
Wider XL armrests	1470086	✓	✓	x	x	✓	✓
Safety lock for XL armrest	1470077	✓	✓	x	x	✓	✓
Safety lock 2.0 (for armrests 2.0)	AP1643760	✓	✓	✓	✓	✓	✓
24" wheels (set)	16379	✓	+	✓	+	x	x
Calf and amputee support right	AP1643573	✓	✓	✓	✓	✓	✓
Calf and amputee support left	AP1643574	✓	✓	✓	✓	✓	✓
Calf support with foot support right	AP1643543	✓	✓	✓	✓	✓	✓
Calf support with foot support left	AP1643572	✓	✓	✓	✓	✓	✓
Lateral support	1535077	✓	✓	✓	✓	✓	✓
Ergonomic headrest holder	AP1531613	✓ <sup>2</sup>	✓ <sup>2</sup>	x	x	✓	✓
Ergonomic headrest	1531465	✓ <sup>2</sup>	✓ <sup>2</sup>	x	x	✓	✓
Ergonomic neckrest	1531608	✓ <sup>2</sup>	✓ <sup>2</sup>	x	x	✓	✓
Toilet pan	10230	✓	✓	✓	✓	✓	✓
Toilet pan holder	AP1602078-XB	✓	✓	✓	✓	✓	✓
Front and back anti-tipper Ergo VIP models	AP1603992	x	x	x	x	✓	✓
Front anti-tipper Ocean Ergo models	AP1601986	✓	✓	✓	✓	x	x
Comfort soft backrest VIP Ergo models	AP1603606	x	x	x	x	✓	✓
Comfort soft backrest Ocean Ergo models	AP1604136	✓	✓	x	x	x	+
Comfort soft backrest Ocean Ergo XL models	AP1604016	x	x	✓	✓	x	x
Tilt-stop screw set	AP60129156	x	x	x	x	✓	✓

1/ only together with protective covers, item no. AP1605002

2/ headrest or neckrest only in combination with Ocean VIP Ergo backrest frame, item no. AP1606985

**New Zealand:**

Freephone 0800 468 222,  
 email [sales@invacare.co.nz](mailto:sales@invacare.co.nz)  
 or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:**

Freephone 1800 460 460,  
 email [orders@invacare.com.au](mailto:orders@invacare.com.au)  
 or visit [www.invacare.com.au](http://www.invacare.com.au)

